

# LUCAS<sup>®</sup> 3 Chest Compression System



## Myth Busted! The LUCAS device does fit large patients!

**Weight:**

180 lbs / 82 kg

**Height:**

6 ft / 183 cm

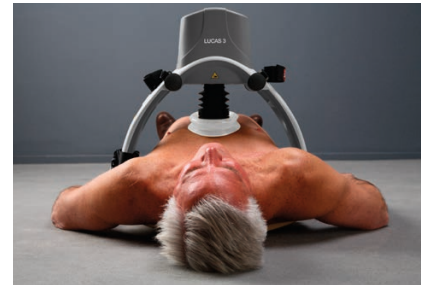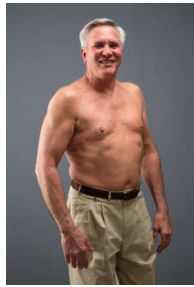


**Weight:**

330 lbs / 150 kg

**Height:**

6 ft, 5 in / 196 cm

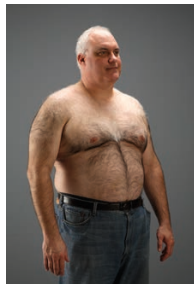


**Weight:**

320 lbs / 145 kg

**Height:**

5 ft, 10 in / 178 cm

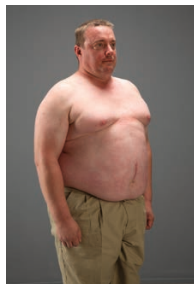

**Proper placement is on upper chest**

Do not assume the LUCAS Chest Compression System can't fit larger patients. Many large patients carry their weight through their mid-section. When LUCAS is properly applied, it rests over the upper chest or sternum, with enough room for the compressor to give compressions according to guidelines. As long as the unit snaps into the back plate and the suction cup does not compress the patient's chest in the start position, it will operate as intended.

**Patients eligible for treatment:**

- Sternum height of 6.7 to 11.9 inches / 170 to 303 mm
- Maximum chest width: 17.7 inches / 45 mm

Using the LUCAS device is not restricted by patient weight.

# LUCAS<sup>®</sup> 3 Chest Compression System

**Myth Busted! The LUCAS device does fit large patients!**



**Weight:**

320 lbs / 145 kg

**Height:**

5 ft, 10 in / 178 cm

**Physio-Control is now part of Stryker.**

**For further information, please contact your local Stryker or Physio-Control representative at 800.442.1142 or visit our website at [www.physio-control.com](http://www.physio-control.com)**

**Physio-Control Headquarters**

11811 Willows Road NE  
Redmond, WA 98052  
[www.physio-control.com](http://www.physio-control.com)

**Customer Support**

P. O. Box 97006  
Redmond, WA 98073  
Toll free 800 442 1142  
Fax 800 426 8049

**Physio-Control Canada**

Physio-Control Canada Sales, Ltd.  
45 Innovation Drive  
Hamilton, ON  
L9H 7L8  
Canada  
Toll free 800 668 8323  
Fax 877 247 7925



Jolife AB, Scheelevägen 17, Ideon Science Park, SE-223 70 LUND, Sweden