



# LUCAS<sup>®</sup> 2|3

Chest Compression System

GENUINE ACCESSORIES FROM PHYSIO-CONTROL

LUCAS 2|3 Chest Compression System





## Thank you for choosing Physio-Control as your partner in helping save lives and improve patient care.

The LUCAS 3 chest compression system is shipped with one battery (dark grey), patient straps, stabilization strap, two suction cups, a hard shell carrying case and the Instructions for Use.

The LUCAS 2 chest compression system is shipped with one battery (light grey), patient straps, stabilization strap, three suction cups, a soft carrying case and the Instructions for Use.

In this catalog you will find the overview of accessories designed to meet your needs and optional power solutions. If you have any questions about accessories or disposables, please contact your Physio-Control representative.

## LUCAS 2|3 Chest Compression System

# Accessories for LUCAS 3 devices

---



### LUCAS Carrying Case, Hard Shell

The LUCAS Hard Shell Case with Top Window for quick battery check and Rear Charge Port is the best way to protect your LUCAS device and carry your spare battery and accessories. Compatible with LUCAS 3 and LUCAS 2 devices.

**11576-000081**



### LUCAS Back Plate, Slim

The slim, low profile back plate is easier to deploy, adjust, strap to transportation devices and clean. Compatible with LUCAS 3 and LUCAS 2 devices.

**11576-000088**



### LUCAS Bumper integrated shaft seal, Black pair

The black rubber bumper accessory matches the colors of the LUCAS 3 device and is available to protect the hinge end caps.

**11576-000091**



### Anti-slip tape, LUCAS Back Plate, Slim

The anti-slip tape is designed to increase the friction of the back plate and to reduce the risk of the device slipping when in contact with other hard surfaces (e.g., spine boards). Six strips in each pack.

**11576-000089**

**11576-000090** (3-Pack)

## Accessories for LUCAS 2 devices



### LUCAS 2 Carrying Bag

LUCAS 2 device is stored and carried in a soft padded carrying bag complete with storage pouches for a spare battery and external power supply.

**11576-000038**



### LUCAS 2 Back Plate

Replacement LUCAS 2 Back Plate is available in the event of damage or contamination of original.

**21996-000044**



### LUCAS 2 Back Plate Anti-slip Tape

The LUCAS 2 Back Plate anti-slip tape is designed to increase the friction of the back plate and to reduce the risk of the device slipping when in contact with other hard surfaces (e.g., spine boards).

**11576-000052**

**11576-000053** (3-Pack)



### LUCAS 2 Rubber Bumper with Integrated Shaft Seal (Pair)

If the hinges protected by the rubber bumper on your device support legs match Fig 1, please order the Rubber Bumper with the integrated shaft seal (pair). The majority of LUCAS 2 devices require this bumper.

**11576-000072**

If you are still unsure which part you need, please contact customer support.




### LUCAS 2 Rubber Bumper (Pair)

If your device has a plastic cap on the support leg hinges matching Fig 2, please order the Rubber Bumper without integrated shaft seal (pair). LUCAS 2 devices older than 2012 may require this bumper.

**11576-000070**

If you are still unsure which part you need, please contact customer support.

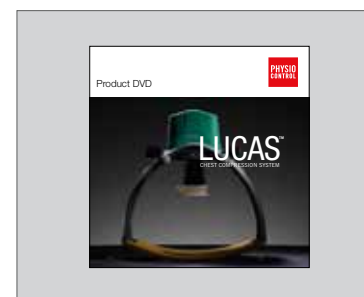

### Product DVD

The LUCAS product DVD features an introduction on the physiology of CPR with Tim Phalen followed by a LUCAS 1 Inservice video and a LUCAS 2 Inservice video. A great DVD to have for new or refresher training.

**21996-000043** (Inservice DVD)

Can also be found on [physio-control.com/LUCAS](http://physio-control.com/LUCAS)

## LUCAS 2|3 Chest Compression System

# Shared Disposables, Batteries & Power Accessories

For use with LUCAS 3 or LUCAS 2 devices



### LUCAS Disposable Suction Cup

LUCAS device has a disposable suction cup which can be easily changed after each use.

**11576-000046** (3-Pack)  
**11576-000047** (12-Pack)



### LUCAS Stabilization Strap

The cushion (new easier to clean material) and straps are used to provide further stabilization relative to the patient and to prevent downward movement of the LUCAS device during operation. Each package includes the neck strap and a pair of support leg straps.

**21576-000074**  
**21576-000075** (4-Pack)



### LUCAS Patient Straps

The LUCAS patient straps are used to secure the patients arms to the support legs of the LUCAS device during transportation.

**11576-000050** (Pair)  
**11576-000051** (Pair, 3-Pack)



### LUCAS Battery

The LUCAS battery is a rechargeable Lithium-Polymer (LiPo) battery which typically last for 45 minutes of operation. You can charge the battery quickly, less than 2 hours in the device if connected to the external power supply or car power cable, and less than 4 hours in the desktop charger.

**11576-000080** (Dark Grey-LUCAS 3 device)  
**11576-000039** (Light Grey-LUCAS 2 device)

Either battery will work in either device, color is the only difference.



### LUCAS Battery Desk-Top Charger

The external battery desk-top charger is a stand-alone unit intended to charge one LUCAS Battery. It typically takes less than 4 hours to charge the battery. A mounting bracket is included for space-saving installation. New dark grey color.

**11576-000060** (U.S., Canada)  
**11576-000061** (EU)  
**11576-000062** (UK)  
**11576-000063** (Australia)  
**11576-000068** (Japan)



### LUCAS Aux Power Supply (S/N 3009 0181 and above)

The external power supply (100-240VAC, 50/60 Hz) can be connected to the LUCAS device and supports prolonged operation as well as charging of the LUCAS Battery while in the device.

**11576-000071** (U.S., Canada)  
**11576-000056** (EU)  
**11576-000057** (UK)  
**11576-000058** (Japan)  
**11576-000059** (Australia)  
**11576-000067** (Switzerland)



### LUCAS Car Cable

The car cable (12-24V DC) can be connected to the LUCAS 2 or LUCAS 3 devices. It supports prolonged operation as well as less than 2 hour charge time of the LUCAS 2 or LUCAS 3 Battery while in the device.

**11576-000048**



### LUCAS PCI Back Plate

The carbon fiber LUCAS PCI Back Plate is intended specifically for use in the cath lab. It is fully radiotranslucent, with minimum shadows. The PCI Back Plate is compatible with LUCAS 2 and LUCAS 3 devices.


**11576-000064**

# Instructions for Use and Training



## LUCAS 3 Instructions for Use - (Software Version 3.0)

26500-003716 English (EN)	26500-003734 Turkish (TR)
26500-003717 International-English (INTL-EN)	26500-003735 Croatian (HR)
26500-003718 German (DE)	26500-003736 Chinese (ZH)
26500-003719 French (FR)	26500-003737 Slovakian (SK)
26500-003720 Dutch (NL)	26500-003738 Korean (KO)
26500-003721 Danish (DA)	26500-003739 Hebrew (HE)
26500-003722 Swedish (SV)	26500-003740 Bulgarian (BG)
26500-003723 Norwegian (NO)	26500-003741 Romanian (RO)
26500-003724 Finnish (FI)	26500-003742 Estonian (ET)
26500-003725 Italian (IT)	26500-003743 Latvian (LV)
26500-003726 Greek (EL)	26500-003744 Lithuanian (LT)
26500-003727 Spanish (ES)	26500-003745 Slovenian (SL)
26500-003728 Portuguese (PT)	26500-003746 Icelandic (IS)
26500-003729 Czech (CS)	26500-003747 Serbian (SR)
26500-003730 Polish (PL)	
26500-003731 Hungarian (HU)	
26500-003732 Russian (RU)	
26500-003733 Japanese (JA)	

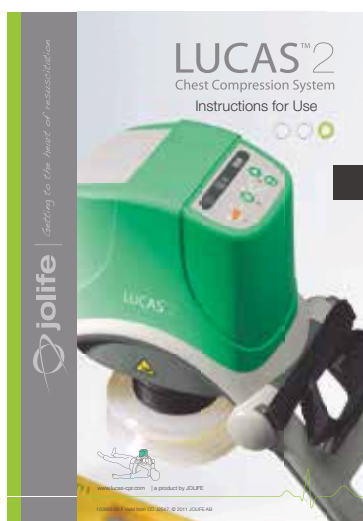

## LUCAS 2 Instructions for Use - (Software Version 2.2)

26500-003588 English (EN)	26500-003606 Turkish (TR)
26500-003589 International-English (INTL-EN)	26500-003607 Croatian (HR)
26500-003590 German (DE)	26500-003608 Chinese (ZH)
26500-003591 French (FR)	26500-003609 Slovakian (SK)
26500-003592 Dutch (NL)	26500-003610 Korean (KO)
26500-003602 Danish (DA)	26500-003648 Hebrew (HE)
26500-003593 Swedish (SV)	26500-003649 Bulgarian (BG)
26500-003594 Norwegian (NO)	26500-003650 Romanian (RO)
26500-003595 Finnish (FI)	26500-003651 Estonian (ET)
26500-003596 Italian (IT)	26500-003652 Latvian (LV)
26500-003597 Greek (EL)	26500-003653 Lithuanian (LT)
26500-003598 Spanish (ES)	26500-003654 Slovenian (SL)
26500-003599 Portuguese (PT)	
26500-003600 Czech (CS)	
26500-003601 Polish (PL)	
26500-003603 Hungarian (HU)	
26500-003604 Russian (RU)	
26500-003605 Japanese (JA)	

For version 2.1 please contact customer service.



## LUCAS 3 or LUCAS 2 Training Device

The LUCAS 3 and LUCAS 2 Training Device is a non-clinical device intended for training purposes only. It should only be used on manikins.

LUCAS 3 99576-000042

LUCAS 2 99576-000020

Lisätietoja on saatavissa paikalliselta Physio-Controlin edustajalta tai internetistä osoitteesta [www.physio-control.fi](http://www.physio-control.fi)



**Physio-Control Headquarters**  
11811 Willows Road NE  
Redmond, WA 98052 USA  
Tel 425 867 4000  
[www.physio-control.com](http://www.physio-control.com)

**Physio-Control Operations  
Netherlands B.V.**  
Galjoenweg 68  
6222 NV Maastricht  
The Netherlands  
Tel +31 (0)43 3620008  
Fax +31 (0)43 3632001

**Physio-Control Finland Sales Oy**  
[info.nordic@physio-control.com](mailto:info.nordic@physio-control.com)  
Tel +358 9 876 46 30  
[www.physio-control.fi](http://www.physio-control.fi)